



Magnifiers

Wave LED

A larger viewing area and LED energy saving technology.

The Wave LED has a new sleek head design to complement any office décor. Offered in a unique and powerful 3.5 diopter magnifier, the Wave LED is superior in both form and function. The dual independent LED light sources offer enhanced performance, as well as a user controlled energy saving option.

Features


- Shadow-free magnification using symmetrical lighting or 3-D magnification using asymmetrical magnification.
- Self-balancing shade and hands-free neck assembly make precise adjustments easy.
- 45" articulating K-arm standard.
- 6.75" x 4.5" rectangular white crown optical-quality glass lens.
- Equipped with two 6W individually controllable LED modules.
- Automatic 9/4 hour auto shut-off feature and (0-50-100%) step dimming to support energy savings.
- No drift spring-balanced arm allows complete freedom of movement.
- Standard table edge clamp mount included. Optional accessories available: surface mount bracket and weighted rolling floor stand.




PHILIPS
burton



Wave LED

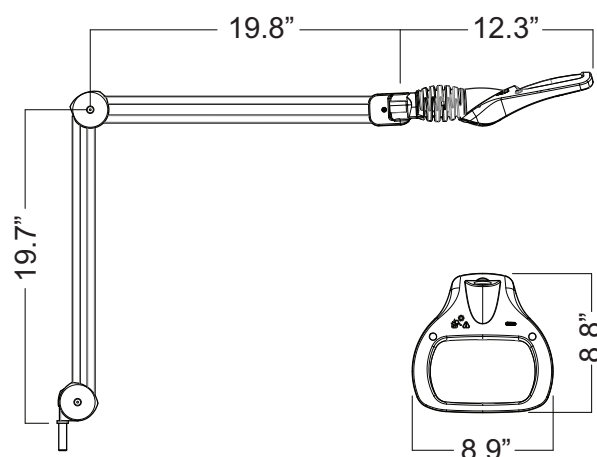

Optional Accessories

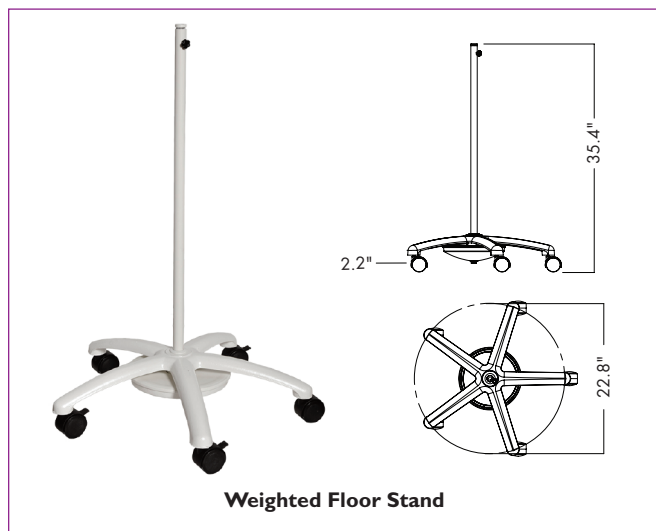

Summary of technical data

Illuminance:	4600 lux (typical)
Color Temperature:	4000 K (typical)
CRI (Color Rendering Index):	80
Light Source:	Two 6W dimmable LED modules
Rated Life of Light Source:	50,000 hours average LED life
Power Consumption:	15 Watts (typical)
Focal Length:	11"
Optics:	3.5 (1.88X) diopter, 6.75" x 4.5" rectangular white crown optical quality glass lens
Dimming:	Step dimming 0 – 50 – 100%
Timer:	9/4 hour auto shut off
Body Material & Color:	Steel arm, die-cast aluminum shade, fully enclosed neck
Arm Technology & Movement:	Heavy duty internal spring 45" parallel three-pivot K-arm
Mounting:	Weighted floor stand, edge clamp or surface mount
Certifications:	UL/cUL listed
Product Warranty:	1-year warranty

Ordering Information

I 20V Model	Description
WML35W45EC120	Wave LED (standard table edge clamp mount included)
Optional Accessories	
WMLWFS	Weighted Floor Stand
WMLSMB	Surface Mount Bracket
WML4DS	4 Diopter STAYS Lens
WML6DS	6 Diopter STAYS Lens
WML10DS	10 Diopter STAYS Lens

**All products available in white only.*



© 2015 Philips Burton. All rights reserved. Product designs protected by copyright. Philips reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

LIT-0091 Rev I



Philips Burton
11500 Melrose Avenue
Franklin Park, IL 60131
Phone: (800) 444-9909
Fax: (847) 451-0744
www.burtonmedical.com